


ペニシリン系抗生物質
Penicillin antibiotics

U250318A

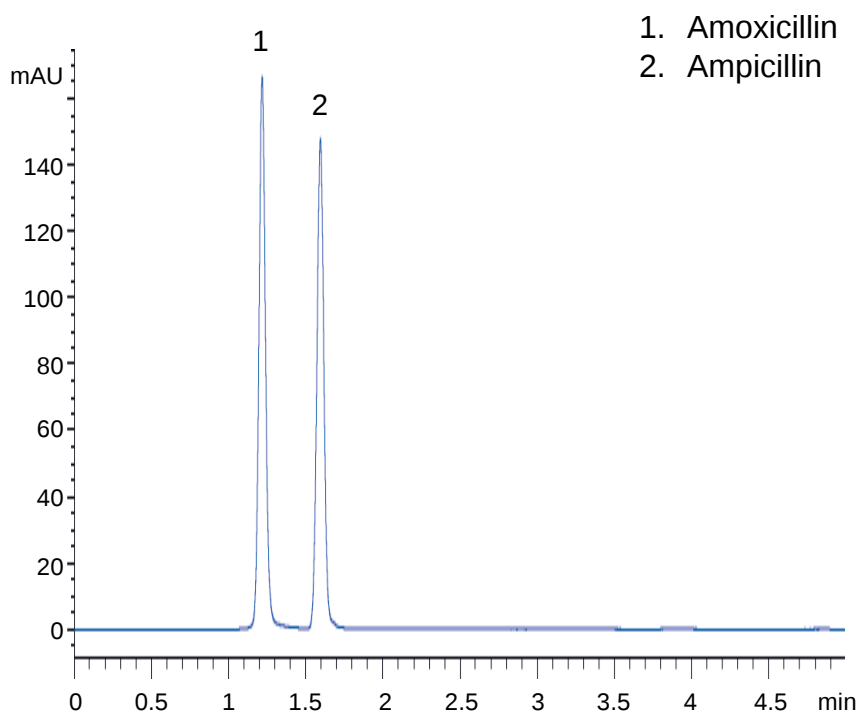

For more details, please see reference. [Rapid-HPLC-analysis-of-Amoxicillin-and-Ampicillin-with-YMC-Triart-C18-ExRS](#)

Column	: YMC-Triart C18 ExRS (5 μ m, 8 nm) 150 X 4.6 mmI.D.
Eluent	: 50 mM sodium phosphate buffer (pH 4)/acetonitrile (82/18)
Flow rate	: 1.0 mL/min
Temperature	: 25°C
Detection	: UV at 220 nm
Injection	: 5 μ L (each 0.25 mg/mL)